



# PuraStat<sup>®</sup>

ABSORBABLE HAEMOSTATIC AGENT

*Pre-filled  
peptide-based  
biomaterial  
for clear  
haemostatic  
control<sup>1</sup>*



PuraStat® is a different kind of haemostat. It is a peptide based haemostatic agent designed to create a 3-dimensional scaffold as a physical barrier that provides atraumatic haemostasis. It is ready without delay. Available in a prefilled ready to use syringe with no preparation required.

### PuraStat® The science behind PuraStat®

PuraStat® is a self-assembling resorbable peptide aqueous solution. The 2.5% aqueous peptide solution is composed of 4 repeating amino acids: R: Arginine, A: Alanine, D: Aspartic Acid -- so the final peptide in PuraStat® is RADA16.

When exposed to physiological fluids such as blood, the aqueous solution neutralises or becomes alkaline. As a result, the  $\beta$  structure peptide molecules self-assemble into nanofibers on a scale similar to the extracellular matrix, and as a result the peptide molecules quickly form a three dimensional meshwork of nanofibers in the aqueous solution, yielding a peptide hydrogel.


When PuraStat® is applied at the point of bleeding, this ionic exchange and pH balance is instantaneous causing a physical (mechanical) closing of the superficial part of the broken blood vessel, resulting in haemostasis (induces natural blood coagulation).

PuraStat® is fully chemically synthesised and does not contain any animal or human derived substance: and any PuraStat® remaining in the body after haemostasis has been achieved, is absorbed no longer than 30 days after application.

### PuraStat® Mode of Action

PuraStat® is a slightly viscous solution of synthetic peptides. Contact between PuraStat® and blood causes the acidic peptide solution to be neutralized and exposed to ions, resulting in the formation of  $\beta$ -sheets

that then form a 3-dimensional scaffold structure that is instantaneous, causing a physical (mechanical) closing of the superficial part of the broken blood vessel, resulting in haemostasis (induces natural blood coagulation).

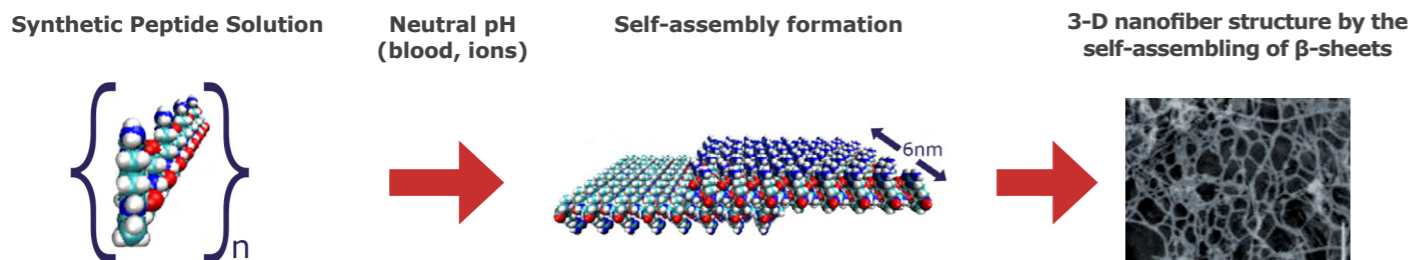


## Procedural Applications

### Prostatectomy Open, Laparoscopic and Robotic

#### Potential Bleeding Sites


- Prostatic fossa along neurovascular bundles
- Vascular pedicles
- Tissue adherent to the tumor
- Tissue surrounding the prostate
- Along the urethrovesical anastomosis

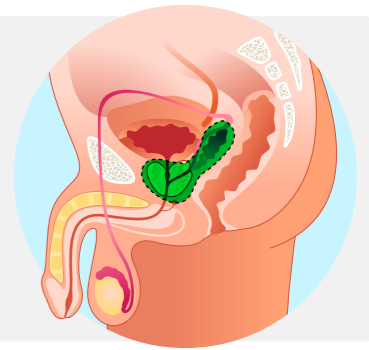

### Partial Nephrectomy

#### Potential Bleeding Sites


- Renal tissue/capsule following tumor removal
- Suture line oozing
- Soft tissue oozing

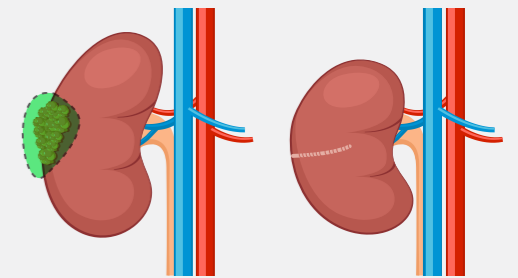

### Cystectomy & Ileal Conduit

#### Potential Bleeding Sites

- Neurovascular bundles
- Seminal vesicles structures
- Vascular pedicles
- Tissue adherent to tumor
- Oozing on suture lines
- Anastomotic suture line ooze

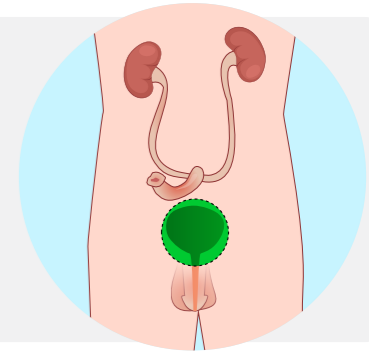

## PuraStat® Advantages




### Transparent Haemostat

- Maintains clear sight of bleeding point



### Reapplicable & Removable

- Can be used repetitively during the same procedure
- Can be removed easily if desired

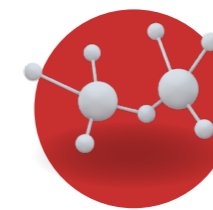

### Synthetic Peptide

- Inert Material
- Well accepted by the body
- Non-swelling<sup>2</sup>



### Easy to Use

- Single prefilled, ready-to-use syringe
- No preparation required
- Does not obstruct the applicator<sup>3</sup>



### Suitability

- Suitable for endoscopic/laparoscopic use
- Can be used in combination with clips, sutures and cautery
- Applicable to narrow spaces
- Covers uneven surfaces

# Illustration of Laparoscopic application of PuraStat®

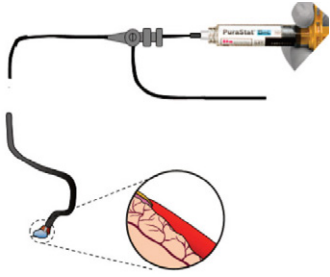


Fig. 1: Oozing bleeding from tissue (parenchyma a/o blood vessels)

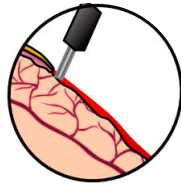


Fig. 2: Remove as much blood as possible from haemorrhagic site (using suction a/o gauze)

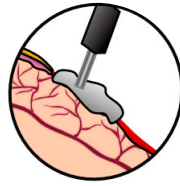


Fig. 3: Apply PuraStat® as close as possible to the bleeding point (also when potentially reapplying)

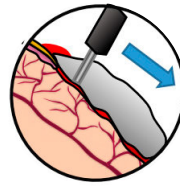


Fig. 4: Continue to apply PuraStat® by moving applicator until the product exceeds the margins of the lesion


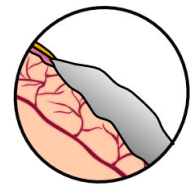


Fig. 5: Haemostasis is achieved

## PuraStat® Nozzle System Type-L (Laparoscopic Catheter)

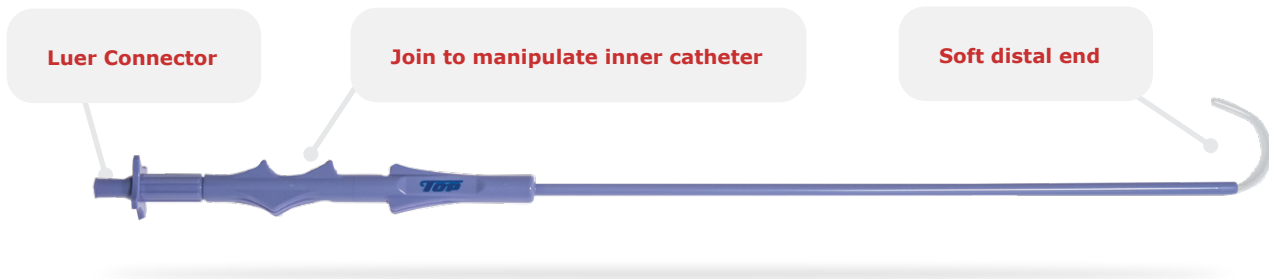

**Consists of 2 layers:** stiffer inner layer and soft outer layer, creating a soft tip

Flexible and curved tip for precise application

When pushed out, the soft tip of the inner tube articulates, enabling accurate application of PuraStat®

## Ordering Information



### Customer Service

australia@3d-matrix.co.jp  
1800 3DMATRIX (1800 336 287)

Product Code	Product Name	Package
621-014	PuraStat® 3mL	1 each
621-015	PuraStat® 5mL	1 each
634163	PuraStat® Nozzle System Type-L	1 box (5 each)
634142	PuraStat® Nozzle System Type-E	1 box (5 each)

### References

1. PuraStat IFU-002 Rev 2.2 2019/04
2. Data on file (Swell Report 2018, Eun Seok Gil 12/06/2018)
3. Risk Management report

### PuraStat®



3-D Matrix Europe SAS  
11 chemin des Petites Brosses  
69300 Caluire-et-Cuire-FRANCE  
Tel: +33 (0)4 27 19 03 40  
infoeu@puramatrix.com  
www.3dmatrix.com



Availability of these products may vary from country or region to another, as a result of specific local regulatory approval or clearance requirements for sales in such country or region. Please contact 3-D Matrix Europe should you need any additional information on device classification.

PuraStat® is not available in the US

### PuraStat® Nozzle System (Catheter)



TOP Corporation  
19-10 Senju-nakai-cho  
Adachi-ku, Tokyo  
120-0035 Japan



Availability of these products may vary from country or region to another, as a result of specific local regulatory approval or clearance requirements for sales in such country or region.

### LOCAL DISTRIBUTOR

Caribbean Park, 44 Lakeview Dr  
Scoresby VIC 3179  
Australia  
e: australia@3d-matrix.co.jp  
p: 1800 3DMATRIX  
(1800 336 287)

